










Reusable Fabric Masks Superfox AM/2S

*Swiss Community Masks, produced according to the requirements of the **Swiss National Covid-19 Science Task Force**.*

-  *Reusable: Guaranteed washable up to 20 times at 60°*
-  *According to the guidelines of NCS-TF*
-  *Skin friendly textiles with OEKO TEX STANDARD 100*
-  *Cup style cut: comfortable up to 8 hours of wearing*
-  *With nose wire: comfortable fit when wearing glasses*
-  *3 different sizes and adjustable earloops*
-  *Water and moisture proof fabrics (ISO 9237, EN 14683: 2019 + AC: 2019)*

**Engineered &
designed in
Switzerland+
by Superfox**

OEKO-TEX®
CONFIDENCE IN TEXTILES
STANDARD 100 



SOLID COLORS

Royal Blue / SFMSK001



Deep Navy / SFMSK002



Violet / SFMSK003



Burgundy / SFMSK004



Mustard / SFMSK005



Burnt Orange / SFMSK006



Pastel Green / SFMSK007



Dark Mint / SFMSK008



Turquoise / SFMSK009



Clear Sky / SFMSK010



Peach / SFMSK011



Magenta / SFMSK012



Light Lime / SFMSK013



Coral / SFMSK014



Snow White / SFMSK015



Black / SFMSK016



SWISS DESIGNS

Swiss Cross / SFMSK017



Swiss Pattern / SFMSK018



Swiss Red / SFMSK019



Swiss Inverse / SFMSK020



Edelweiss Red / SFMSK021



Edelweiss White / SFMSK022



Paisley Red / SFMSK023



Paisley White / SFMSK024



VARIOUS MODELS

Roses / SFMSK025



Roses Baroque / SFMSK026



Hibiscus / SFMSK027



Ice Rose / SFMSK028



Lemons / SFMSK029



Tangerine / SFMSK030



Lemon Slice / SFMSK031



Mojito / SFMSK032



Cats / SFMSK033



Dogs / SFMSK034



Fox / SFMSK035



Jellyfish / SFMSK036



Jungle / SFMSK037



Tiger / SFMSK038



Forest Animals / SFMSK039



Superfox / SFMSK040



Pink Unicorn / SFMSK041



Magic Unicorn / SFMSK042



Blue Sky / SFMSK043



Mermaid / SFMSK044



Space Animals / SFMSK045



Monsters / SFMSK046



Boho Moon / SFMSK047



Succulents / SFMSK048



Camouflage 1 / SFMSK049



Camouflage 2 / SFMSK050



Camouflage 3 / SFMSK051



Camouflage 4 / SFMSK052



Optical Illusion / SFMSK053



Colors Geometry / SFMSK054



Leaves B/W / SFMSK055



Colors Gradient / SFMSK056



Japan / SFMSK057



Rainbow / SFMSK058



Ethno / SFMSK059



Ice-cream / SFMSK060



Astronaut / SFMSK061



Vintage Astronaut / SFMSK062



Adventure / SFMSK063



Mountain / SFMSK064



Music / SFMSK065



Notes / SFMSK066



Origami / SFMSK067



Surf / SFMSK068



Skulls Black / SFMSK069



Skulls Gold / SFMSK070



Halloween / SFMSK071



Tropical Abstract / SFMSK072



Autumn Trees / SFMSK073



Japanese Waves / SFMSK074



Matisse Dance / SFMSK075



Watercolor View / SFMSK076



Smiley / SFMSK077



Snake / SFMSK078



Mistletoe / SFMSK079



Winter Adventure / SFMSK080



SEASONS - AUTUMN

Bushes / SFMSK081



Acorns / SFMSK082



Mushrooms / SFMSK083



Hedgehog / SFMSK084



Autumn Leaves / SFMSK085



Sox / SFMSK086



Bouquet / SFMSK087



Acorn / SFMSK088



SEASONS - WINTER

Winter Knit 1 / SFMSK089



Winter Knit 2 / SFMSK090



Winter Knit 3 / SFMSK091



Winter Knit 4 / SFMSK092



Winter Knit 5 / SFMSK093



Ice Skating / SFMSK094



Hats / SFMSK095



Didi / SFMSK096



Winter Forest / SFMSK097



Didi & Dodo / SFMSK098



Ecru Knit / SFMSK099



Orange Knit / SFMSK100



Black Knit / SFMSK101



Mixed Knit 1 / SFMSK102



Mixed Knit 2 / SFMSK103



Mixed Knit 3 / SFMSK104



CHARITY "HEIDI HILFT" FOUNDATION *

Heidi 1 / SFMSK105



Heidi 2 / SFMSK106



* For every mask sold from these series, Superfox.ch donates **CHF 2.** - to the Swiss non-profit organization "**Heidi Hilft**". So your purchase not only promotes health, but also help, faith and hope.

CURLING

I Throw Rocks / SFMSK107



Curling ♡ / SFMSK108



EARLOOP COLORS

White



Black



Light Grey



Dark Grey



Light Blue



Navy Blue



Beige



Rosy Pink



Yellow



Orange



Jungle Green





Gold





SIZES



FACE MASK SIZE L

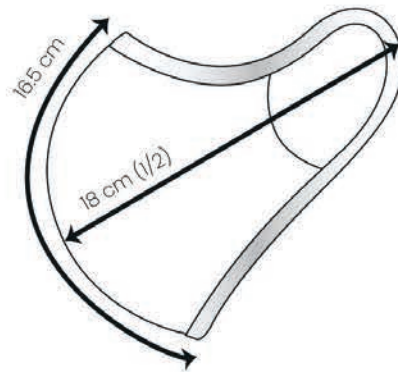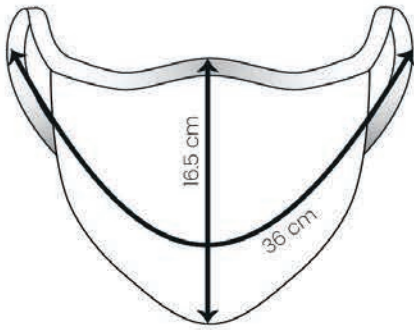


FACE MASK SIZE M

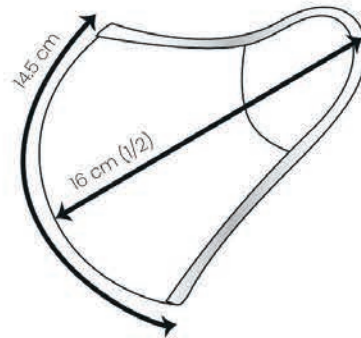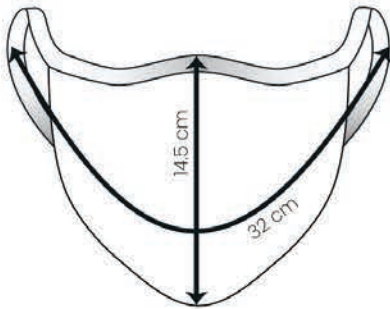

FACE MASK SIZE S

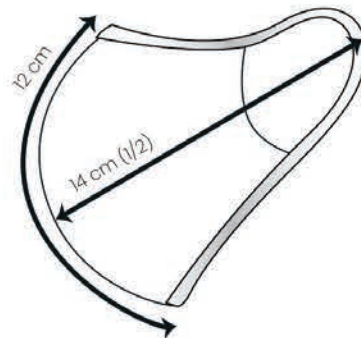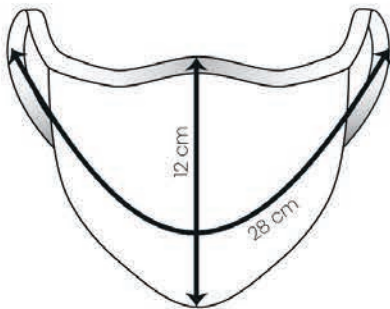

CUSTOM

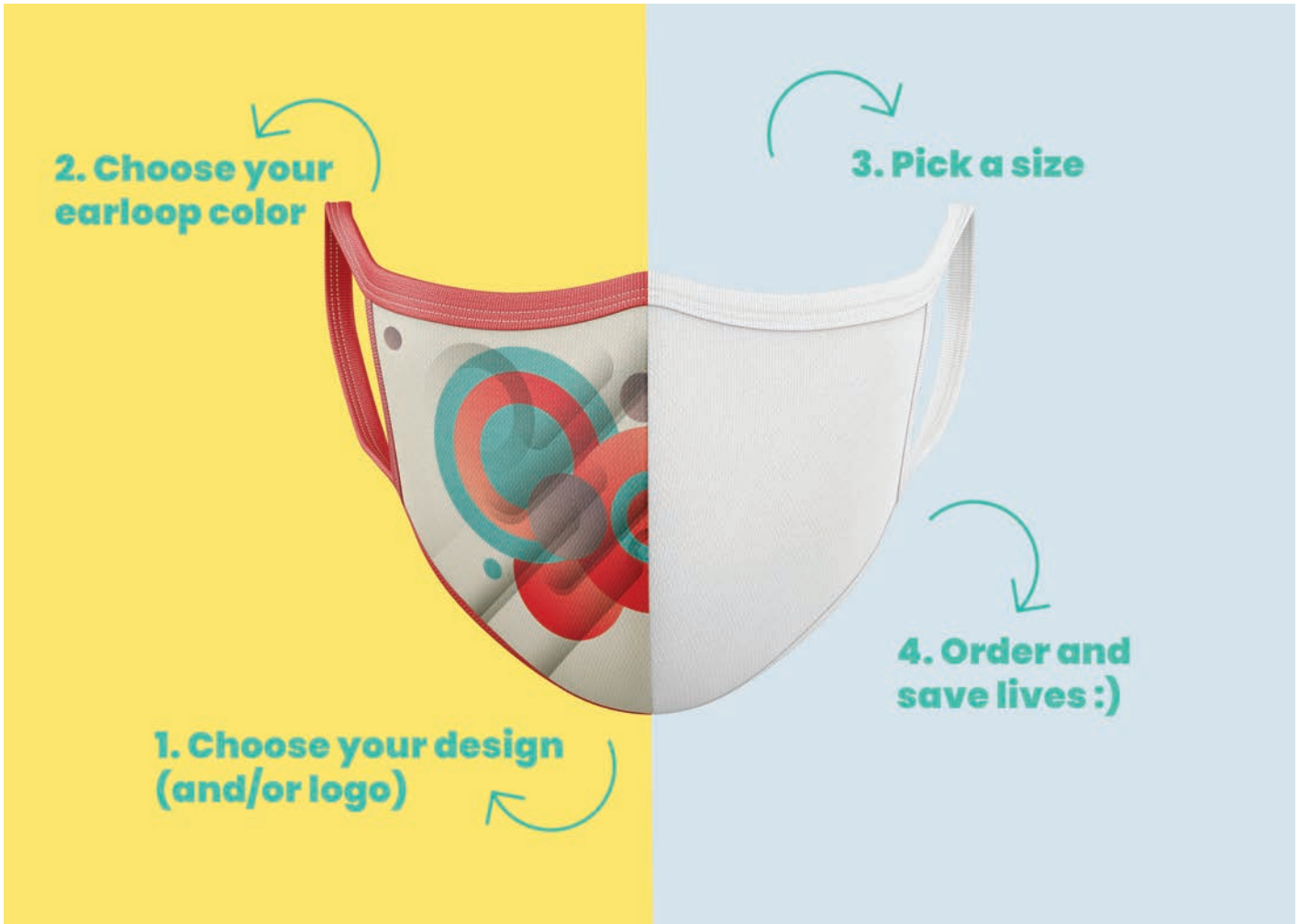

LABEL

REUSABLE RESPIRATOR/HYGIENE FACE MASK MOD. SUPERFOX AM/2S



- Protects airways from aerosoles and moisture drops
- Durable to hard particles
- Waterproof surface
- Respirator: Cup-Style (fully covers nose and mouth)
- Suitable for most face shapes and sizes
- Keeps its shape while wearing thank to inserted nose clip
- Washable with 60°C (tested up to 15 times)
- Can be sterilised in autoclave



COMPLIANCE WITH THE TECHNICAL REQUIREMENTS ACCORDING TO BOEN I09 PROVISION 4525.

Order SND / 354/2020, of April 19, which establishes exceptional measures to guarantee access to the population of recommended use as hygienic measures for the prevention of infections by COVID-19:

It complies with the technical specifications UNE 0064-1: 2020, UNE 0064-2: 2020 or UNE 0065: 2020 that establish the minimum requirements that hygienic, non-reusable and reusable masks must meet in terms of the materials used for their design, elaboration, preparation, marking and use or any equivalent standard.

It also complies with the EN 14683: 2019 + AC: 2019 standard (points 5.2.2, 5.2.3 and 5.2.4) that establishes, according to Directive 93/42 / EEC, the operating requirements for Single-Use Surgical Masks Type I.

Test based on EN 14683: 2019 + AC: 2019
Test report obtained in Aitex accredited laboratory, protocol number 2020TM1188 dated 05/05/2020

Number of test replicates: 5

TESTS	RESULTS /Average ± SD/
Pt 5.2.2. Bacterial Filtration Efficiency (BFE) %	95,06 ± 3,24
Pt 5.2.3. Breathability: Differential Pressure (Pa/cm2)	17 ± 0
Pt 5.2.4. Splash resistance pressure (kPa)	Failure 0 of 3 at 17 kPa

HOW TO USE


If possible wash the mask before using for the first time.

1. Try the mask on and check if the default length of the elastic strings is ok. If the mask is loose adjust the length in with the rubber bands by creating small tight knots.

2. The mask is equipped with a metal clip that has to be upward (placed on your nose).The bottom part of the mask needs to cover the chin. Apply the hooks to your ears.

3. Once the mask is applied, adjust the metal clip to the form of the nose and make sure it is tightly fixed.

4. Stretch/cover the mask with your palms, inhale and exhale. If you feel that there is air coming out from the edges of the mask, readjust the length of the hook according to the pictures.

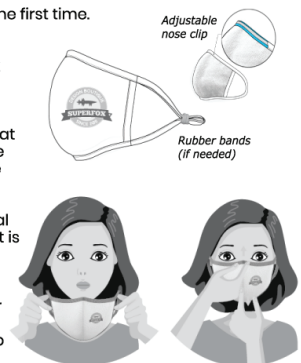

The mask needs to be either washed or sanitized in autoclave after each use - contact with people or use in an environment with aerosols. You must dispose/stop using your mask if its structure is corrupted, affected in any way. Throw away the mask if it can no longer stay tight on your face and air gets out from edges. If you use the mask for more than 8 hours per day you must wash or sterilize it before the next use. Store in a dry and ventilated place.

COMPOSITION

Outside: 100% Polyester | Filters: 100% Polypropylen
Tested and certified with EN ISO 4920: 2012, EN ISO 9237: 1999, EN ISO 6330: 2012

Made in Bulgaria (EU) for Univex AG, Haus zum Hirschen, 8758 Obstalden, Switzerland

